

Root® with SedLine® Brain Function Monitoring

A More Complete Picture Starts With More Complete Data




Root with SedLine Brain Function Monitoring helps clinicians monitor the state of the brain under anesthesia with bilateral data acquisition and processing of Electroencephalogram (EEG) signals

- > Four simultaneous channels of frontal EEG waveforms
- > Patient State Index (PSI) is a processed EEG parameter that is related to the effect of anesthetic agents
- > A Density Spectral Array (DSA) display, which contains left and right spectrograms representing the power of the EEG on both sides of the brain
- > Multiple screen views, which expand information while enabling customization in the operating room and intensive care unit



Monitor Display

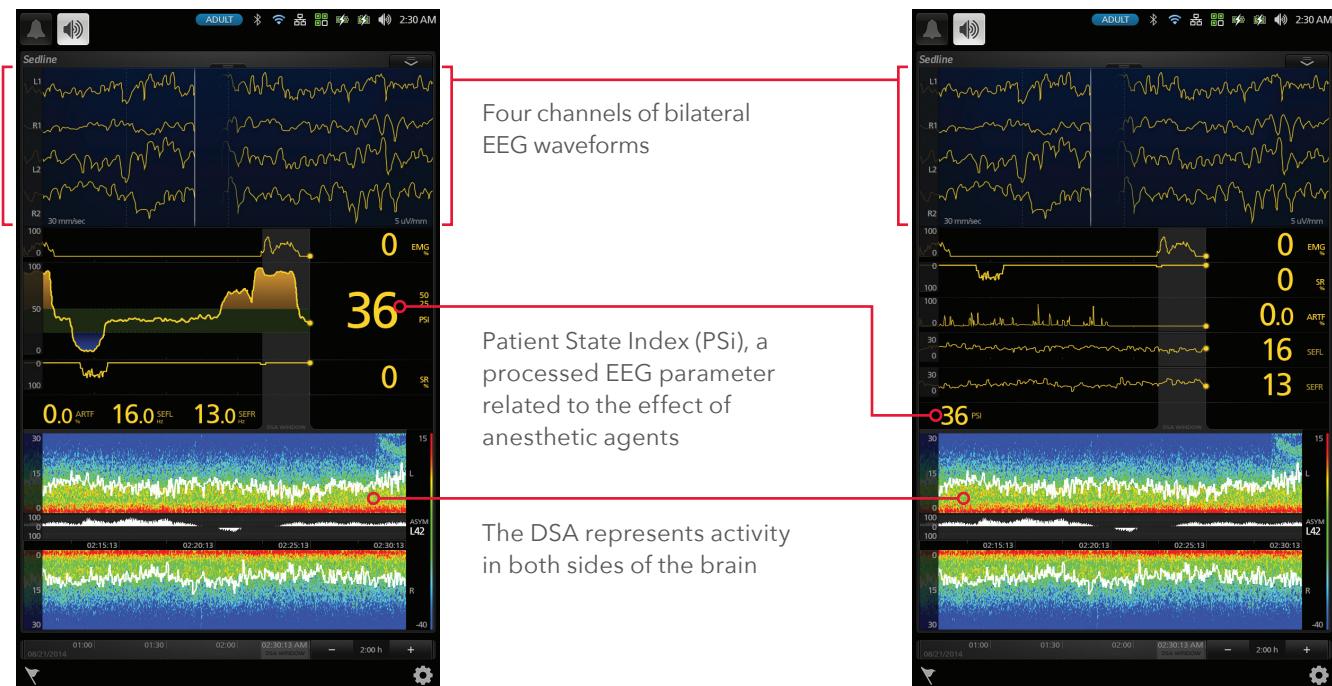

The SedLine EEG Sensor

- > Four active leads collect data from the frontal lobe
- > Allows simultaneous application of SedLine and O3® Regional Oximetry sensors
- > Soft foam pads improve patient comfort



The SedLine module plugs into the Root patient monitoring platform via Masimo Open Connect™ (MOC-9™) ports

SedLine Specifications

PHYSICAL CHARACTERISTICS

Module Physical Dimensions

Width	13/10 in (33mm)
Length	4 in (102mm)
Thickness3/4 in (19mm)

ENVIRONMENTAL

Module Operating Conditions

Temperature at Ambient Humidity 5-40°C

Module Storage and Shipping Conditions

Temperature at Ambient Humidity -40-70°C
Storage Humidity 15-95%, non-condensing
Exposure to Pressure 500-1060 mbar

Sensor Specifications

Patient Weight	> 30 kg
Application Site	Forehead
Active Channels	4
Active Electrodes	L1, L2, R1, and R2

Ground Electrode	CB
Reference Electrode	CT
Duration of Use	Maximum of 24 hours
Latex Content	Does not contain natural rubber latex

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

Masimo U.S.
Tel: 1 877 4 Masimo
info-america@masimo.com

Masimo International
Tel: +41 32 720 1111
info-international@masimo.com

