
Artromick AC Medication Carts

Secure Medication Organization



Artromick Understands Mobile Computing

Artromick brings its superior design expertise to the hospital market with the Artromick AC Medication Cart Series. This innovative product line features the latest design and serves as a platform for existing and emerging medication control technology.

Artromick has redefined medication cart security with a broad choice of lock systems including keyless entry with automatic relocking, mag stripe card access, keyless narcotic storage and cart user audit.

Best of all, Artromick AC Medication Carts offer huge capacity and intelligent solutions for medication management and organization.

We welcome the opportunity to earn your business and look forward to the time when we can illustrate the advantages of our Artromick AC Medication Carts.

AC Medication Cart Highlights

- A Internal pass-thru waste receptacle ensures a trim-line footprint
- B Expansive work area with slide-out writing surface
- C Magnetic stripe, proximity or bar code security for cart access
- D Optional keyless access with auto-relocking
- E Integrated and flexible drawer divider system
- F Configurable medication storage with 2 bin sizes and 3 drawer sizes
- G Full selection of convenience accessories and options
- H Stable, trim-line footprint for simple maneuverability



Cassette & Bin Options

Artromick gives you a choice of 2 or 3 tier durable cassettes with the option of 5.5" or 8" bins.

2-Tier & 3-Tier Cassettes



The AC Medication Cart accommodates 2 & 3 tier cassettes to customize the medication cart configuration that best suits your requirements.

8" & 5" Bins



AC cassettes give you a option of 8" and 5" bins that offer a flexible divider system. Large bin lens accommodates clear patient labeling.

Security Options

Artromick AC Medication Carts offer you a choice of lock systems that best meets the specific needs of the care environment. Our advanced AutoLock™ system offers unparalleled security in a proven lock system designed specifically for healthcare carts. These technology upgrades make this lock system the most advanced cart security available.



AutoLock™

Automatically secures a motionless unlocked medication cart with a user-defined relock time delay from 1-99 minutes. The system also permits the use of over 400 individual user codes to further enhance cart access tracking and security. See available upgrades below.

Keyless Access with Manual Relocking

Simply enter a 4-digit access code, push down the lock lever and the cart is ready for access. To secure the medication cart, simply push in the lock lever. Clear visual indication of lock status is apparent from a distance.

Key Access with Manual Relocking

Featuring the Best® Core Removable Lock System, this basic yet proven lock system provides your medication cart advanced security and the flexibility to change keys simply and easily.

Security Upgrades




Keyless Narcotics Storage

The optional keyless narcotics storage feature greatly augments security to further optimize control and accountability.



Secondary Access Security

Select from a bar code, mag stripe or proximity card reader to serve as the primary or secondary access requirement for the AC Medication Cart.

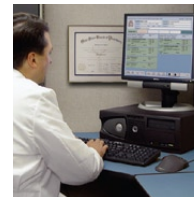

Remote Management

Wirelessly enables authorized users to easily program all carts and retrieve event logs, add / delete user profiles and adjust cart settings from a single PC.




Audit Tracking

Authorized users can track cart activity, including user access and relocking for over 2,000 events.

Accessories & Options


Streamlined nurse workflow is assured with a full complement of convenience accessories, organizing the medication administration process at the point-of-care.

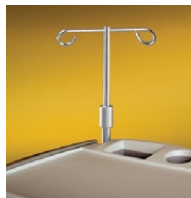

IV Pole




Drawer Dividers

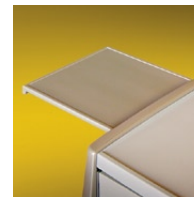

Writing Surface




Sharps / Gloves

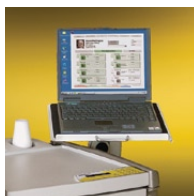

Mobile Computing



Cart Light

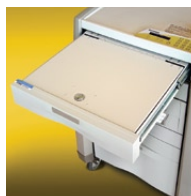


Narcotics Drawer



Waste Container



Storage Module

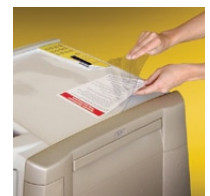

Clear Top Mat

AC Medication Cart Features & Specifications

The Artromick AC Medication Cart line combines the latest medication cart design with plenty of capacity and simple maneuverability. Choose the drawer configuration and cassette size that best meets your requirements for highly-secure and organized medication control.

The Artromick AC is perfect for point-of-care medication administration and as a deliver workstation for your automatic dispensing cabinet.

Features

- Power status display on LCD
- Large retractable writing surface
- Internal pass-thru waste receptacle
- Modular bin / storage drawers
- Expansive storage drawers
- Slot-Lock™ security system secures all partially closed drawers and bins left open up to 1 inch
- Keyless access and auto-relocking with AutoLock option
- Wireless remote security control (optional)
- 5" double ball bearing casters (17.7cm)

Benefits

- Improved workflow and accuracy
- Safe and ergonomic design
- Advanced medication security and control
- Lightweight design for easy maneuverability
- Customized medication / supplies storage

AutoLock Charging System

Input Voltage:	115 VAC (range 90–130 VAC) 230 VAC (range 180–260)
Charger:	Multi-Stage Charger
Battery Type:	12V sealed lead acid with 18 amp hours for extended run times
Short-circuit:	15 amp DC breaker
Status Indicator:	Battery status on LCD display
Approvals:	FCC Class A, Part 15 Approval

Standard Cart Specifications

Base Dimensions:	ACS:	24" x 26"	61cm x 66cm
	ACM:	24" x 38"	61cm x 96.5cm
Work Surface Area:	ACS:	16" x 18"	40.5cm x 46cm
	ACM:	16" x 28"	40.5cm x 71cm
Work Surface Height:	Compact:	39.5"	100cm
	Standard:	43"	109cm
Construction:	Formed steel, aluminum, heavy-duty polymer work surface / end panels		
Finish:	High-density polymer		
Wheels:	Swivel casters with 1 locking; tracking caster optional		

AC Cart Models

Available in 2 sizes and 2 heights:

